

# Aesculap Neurosurgery

## Kerrison Noir®

Detachable Bone Punches



**No** irritating reflections

# Kerrison *Noir*<sup>®</sup> No irritating reflections

## *Black*

---

The Noir<sup>®</sup> surface coating is not only a visual feast, it also prevents effectively any irritating light reflections.

## *Hard*

---

The Noir<sup>®</sup> coating increases the surface hardness, to improve the cutting characteristics and extent the service life.

## *Smooth*

---


The smooth Noir<sup>®</sup> surface ensures optimum gliding properties and enhances the sticking resistance.

## *Detachable*

---

The Kerrison Noir<sup>®</sup> bone punch can be taken apart easily, without any tools – for efficient reprocessing.



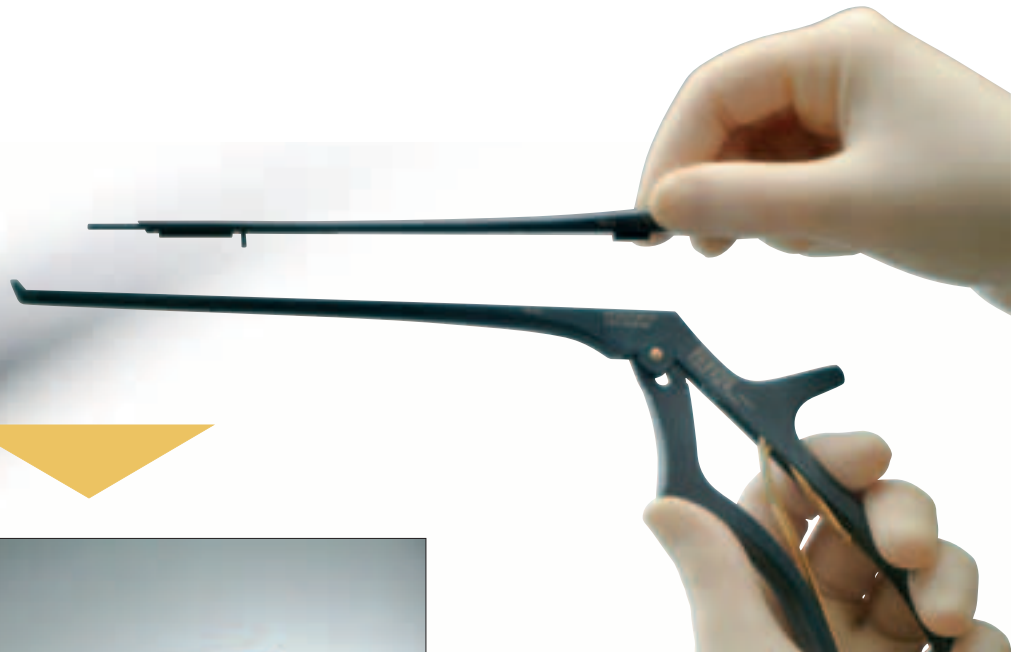

### Easy disassembly


Unhook the spring from the outer handle.



Close the handle completely and pull down the inner handle.



Pull off the slider towards the rear end.

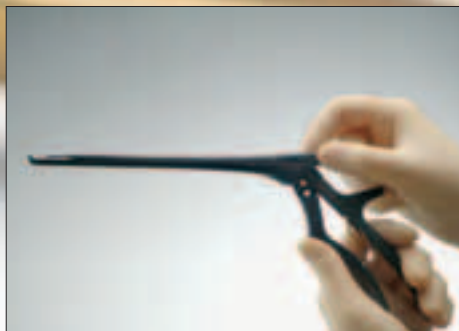

Insert the spring in the outer handle.

Push upwards the inner handle and open the handle.

Attach the slider and push it towards the front end.

### Easy assembly

# Product range

130°  
up



Art. no.	Length	Width	Ejector
FK900B	180 mm, 7"	1 mm	-
FK901B	180 mm, 7"	2 mm	-
FK902B	180 mm, 7"	3 mm	✓
FK903B	180 mm, 7"	4 mm	✓
FK904B	180 mm, 7"	5 mm	✓
FK905B	180 mm, 7"	6 mm	✓
FK913B	200 mm, 8"	2 mm	-
FK914B	200 mm, 8"	3 mm	✓
FK915B	200 mm, 8"	4 mm	✓
FK916B	200 mm, 8"	5 mm	✓
FK917B	200 mm, 8"	6 mm	✓

## Punches with thin foot plate

FK906B	180 mm, 7"	1 mm	-
FK907B	180 mm, 7"	2 mm	-
FK908B	180 mm, 7"	3 mm	✓
FK909B	180 mm, 7"	4 mm	✓
FK910B	180 mm, 7"	5 mm	✓
FK961B	200 mm, 8"	1 mm	-
FK962B	200 mm, 8"	2 mm	-
FK963B	200 mm, 8"	3 mm	✓
FK964B	200 mm, 8"	4 mm	✓
FK965B	200 mm, 8"	5 mm	✓
FK925B	230 mm, 9"	2 mm	-
FK926B	230 mm, 9"	3 mm	✓
FK927B	230 mm, 9"	4 mm	✓
FK928B	230 mm, 9"	5 mm	✓

90°  
up




Art. no.	Length	Width	Ejector
FK948B	180 mm, 7"	1 mm	-
FK949B	180 mm, 7"	2 mm	-
FK950B	180 mm, 7"	3 mm	✓
FK951B	180 mm, 7"	4 mm	✓
FK952B	180 mm, 7"	5 mm	✓
FK953B	180 mm, 7"	6 mm	✓



The Kerrison Noir® punches with a thin foot plate are to be used exclusively for removing tissues and soft, small bone fragments.

They are identified by their golden inner handle.

Kerrison Noir® bone punches with standard footplates have a black inner handle.

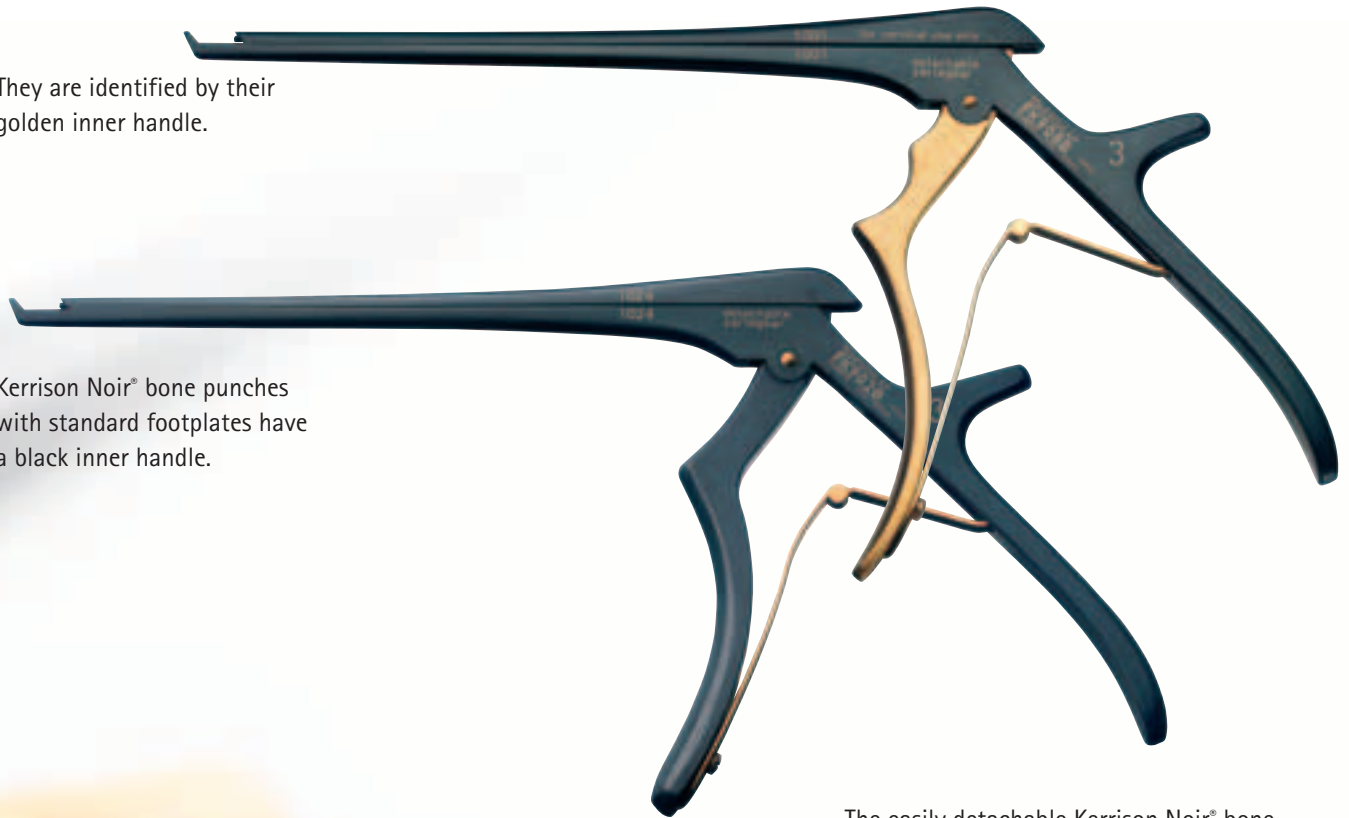

The easily detachable Kerrison Noir® bone punches from Aesculap are identified by their golden spring.

Large number for identification of jaw width at a glance.

Ejector – for the easy removal of punched-out material.

Numerical code – for reliable identification when assembling the two punch components.





# AESCULAP®

## Expertise in Surgery.

**B | BRAUN**  
SHARING EXPERTISE

**Aesculap AG & Co. KG**

Am Aesculap-Platz  
78532 Tuttlingen

Phone +49 7461 95-0

Fax +49 7461 95-2600

[www.aesculap.de](http://www.aesculap.de)

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogues and pricelists and to take legal actions.